

Medi-Pulse Q

SUPERIOR INNOVATION



Combining FOUR modalities in ONE housing

The Medi-Pulse Q offers practitioners the opportunity to dramatically increase their treatment portfolio with no need to purchase multiple devices.

Two highly effective IPL hand pieces offers various treatment options that will provide unrivalled stable power output and efficiency.

A large touch screen interface with built-in pre-set parameters with optimised controls is designed with the practitioner in mind.

With a Nd Yag Q Switch handle for up to 20 treatment options with the 532nm 1064nm and 1320nm spot sizes, the Medi-Pulse Q offers a large range of superior treatments.

Perfect to treat

- Skin rejuvenation and anti-aging
- Skin and pore tightening
- Tattoo removal
- Benign pigmented lesions, deep and superficial
- Age Spots
- Freckles
- Acne scarring
- Active acne
- Carbon peel
- Nail fungus

Skin Rejuvenation
Super Fast Hair Removal
Carbon Laser Facial
Tattoo Removal



Handpieces

Our hand pieces have pure sapphire tips NOT glass. Pure sapphire has a much higher light transmission than glass. For example if you set the energy to 10J/cm you will get 10J/cm, whereas with glass you are more likely to actually get 5J/cm with the same settings. Our Filters also are Pure Sapphire.

Filters

You want to be able to get the most from your IPL SHR device. We give you a choice between different wavelengths to target different skin conditions, as well as hair removal. For example a 430nm filter is great for treating acne, whereas a 690nm filter is perfect for permanent hair reduction.

We also have the 690nm built into our actual HR hand piece where there is no need to apply a filter.

Filters allow you greater flexibility with treating different skin types, so for the HR (Hair Removal SHR) hand piece we also have a 750 filter for quite darker skin tones.


Having this built-in filter in your hand piece guarantees no energy leakage.

For the SR (Skin Rejuvenation) hand piece however we do need different filters to tackle different diseases which use different wavelengths to break down melanin in the skin, eliminating pigmentation, discolouration, smoothing out skin tone, redness, broken capillaries, also calming and repairing acneic skin.

* Please make sure when choosing your device to check your warranty! Not only on your device but also your hand pieces. Make sure the company you choose offers at least 300,000 shot warranty... we do...

Our hand pieces can reach up to 1 million shots.



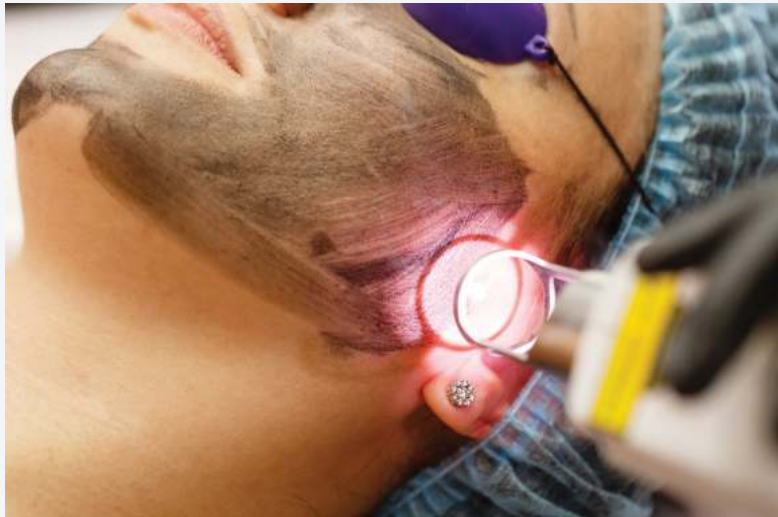

Carbon Peel / Hollywood Facial

The carbon laser facial is a painless, non-invasive treatment. It stimulates collagen production, significantly combats active acne, acne scars, large pores, pigmentation, age spots, melasma and evens out skin tone.

The laser takes a 15 micron layer of skin off, leaving a nice glow with a vibrant complexion. People in Hollywood do this just before the red carpet. No downtime, no discomfort.

AMAZING RESULTS!

There's no such thing as regrets

The Medi-Pulse Q is armed with a Q-switched Nd: YAG 532nm, 1320nm and 1064nm laser for the removal of multi-coloured, professional or amateur tattoos as well as many skin treatments, enabling practitioners to capitalise on the most in-demand treatments.

Operating with intuitive software which allows full control of all parameters, and with easy-to-use detachable laser spot sizes, the Medi-Pulse Q is the perfect choice for all new users of tattoo removal technology.

The Medi-Pulse Q is a robust workstation that will continually deliver results, time after time.

Wavelengths explained

532 Wavelength

This spot size is ideal for treatment of red/orange/coffee ink pigments and superficial pigmentation.



1320 Wavelength

Carbon Laser Peeling for many skin concerns including significant reduction of pore tightening, acne and anti-aging.



1064 Wavelength

For treatments of black ink pigments and deep dermal pigmentation.



Frequently Asked Questions

Tattoo Removal

How many treatments are necessary?

This depends on several factors, including the density and colour of inks and the age of the tattoo. Amateur tattoos can be quick to remove, taking three to eight treatments, whereas professional tattoos vary considerably and some can take up to 10 treatments. Treatment sessions are spaced six to 10 weeks apart.

How does the skin look after treatment?

An instant whitening of the skin is seen after treatment, but this generally fades quickly. The area is likely to feel tender for a few hours and can be red and swollen for up to 48 hours.

Are there any side effects?

Very rarely the skin may become lighter or darker, or a small blister may form depending on sensitivity.

How does the treatment feel?

Most patients describe the sensation as similar to a hot pin prick or the flick of an elastic band. Whilst mildly uncomfortable, most people tolerate the treatment well.

Carbon Laser

How long does the treatment take?

You can easily fit a carbon laser facial into your workday. The entire treatment, including a consultation, the application of the carbon paste and the laser procedure, takes approximately 45 minutes.

What will I look like after treatment?

You could experience some mild redness after the treatment, but you will also see results in your glowing, radiant skin and reduced acne and pore size immediately following your first treatment. This treatment has the power of multiple chemical peels, in a safe, non-invasive laser procedure.

How many sessions will I need?

You will begin to see powerful results after only one session. A course of four to five sessions, spaced one to two weeks apart, can give great results for improved collagen production, reduced sebum and tighter, more radiant skin.

Treatment photos



Above:
Pigmentation Removal

Right:
Super Fast Hair Removal
Results after one session
(4 week interval)

Below:
Carbon Peel treatments



“Quality in a product is not only what you put into it.
It is what the customer gets out of it.”

~ Soma Lasers Australia

Unmatched, incomparable device made with the therapist in mind

Soma devices are made with the highest quality engineering working together with Germany, Korea, Japan the USA and proudly designed in Australia for fantastic durability.

The new technology in this device is seriously unsurpassed. Filter changing is a thing of the past! Everything is programmed for you. In-built filters means more stable power output. Your treatments will blow you away.

Strong, smart, beautiful, a device you can count on for years to come.

Warranty

We back up our promise at Soma Lasers Australia, with a 3 year warranty on all our new medical aesthetic devices.

Marketing

This is a BIG one! What's the point of investing in new technology for your salon or clinic, if you're not attracting any clients? Some of the biggest concerns we hear is 1) How can I feel certain that I will make my money back and 2) I don't know what to post on social media. This is why Soma Lasers Australia INT sets you up with a great kick-start on your marketing! \$2000 Free Marketing Campaign! We guarantee you customers!

YES! We offer you a complete package because we know how daunting it is to set up a successful business!!!



Specifications

Spectrum	HR 640, 690, 755-1200nm SR 430, 480, 530, 590, 640, 690, 755-1200nm
IPL SHR Energy	1-50J/cm ²
Pulse Width	2-10ms
Pulse Delay	5-60ms
SHR Frequency	1-15Hz
Spot Size	HR + SHR 16x50mm ² SR + SSR 10x40mm ²
Flash	1am Xenon Lamp 500,000 max with 300,000 guaranteed
Laser Energy	100-2000mJ
Laser Frequency	1-10Hz
Wavelength	1064/532nm
Voltage	AC100-120V or 220-240V
Dimensions (LxWxH)	500mm x 450mm x 1200mm
Weight	37kg

Equip your salon
or clinic with

soma
Lasers Australia

Our demonstration suite is situated at:
3/10 Atherton Road, Oakleigh Vic 3166

M: +61 403 744 648

E: support@somalasersaus.com.au

W: www.somalasersaus.com.au